


EXHIBIT 73



*Secure to
 conflict:*

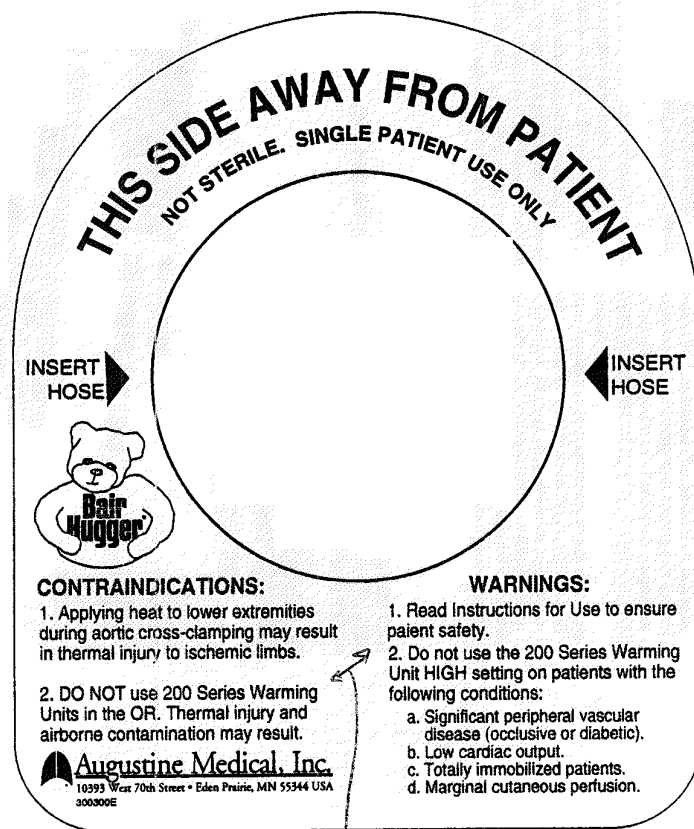

Instructions for Use: Bair Hugger Full Body Surgical Blanket Model 600 Use with any Bair Hugger Warming Unit or the Polar Bair™ Model 600 in the Operating Room.



<p>1</p> <p>Place the Bair Hugger® Blanket on the patient with the perforated side, that is, with small holes, against the patient's skin.</p>	<p>2</p> <p>Insert the Bair Hugger hose into the cardboard Inlet on the Blanket.</p>	<p>3</p> <p>Remove paper backing from the tape barrier and adhere to patient's upper chest/shoulder area. This prevents air migration towards the surgical site.</p>
---	---	---

<p>4</p> <p>Turn the Unit ON and select the appropriate temperature setting. Use this Blanket with any Bair Hugger Warming Unit or the Polar Bair Model 600.</p>	<p>5</p> <p>Place a cotton blanket on top of the Blanket for maximum effectiveness.</p>	<p>6</p> <p>Monitor the patient's temperature on a regular basis, at least every 10 to 20 minutes and adjust the temperature setting of the Unit as required.</p>
---	--	--



Printed on 100% recycled paper, 15% post-consumer. Made in USA.
U.S. Patent Number 4,572,188. Other patents pending © 1994 Augustine Medical, Inc. All rights reserved. Date of issue: 12.1.94

Augustine Medical, Inc.
10000 Augustine Medical Drive
Augustine, CA 95021
408.941.4864 • 800.733.7775 • Fax 408.941.4898